


Service Manual

Domestic Beauty Machine
Warm & Cool Beauty Device



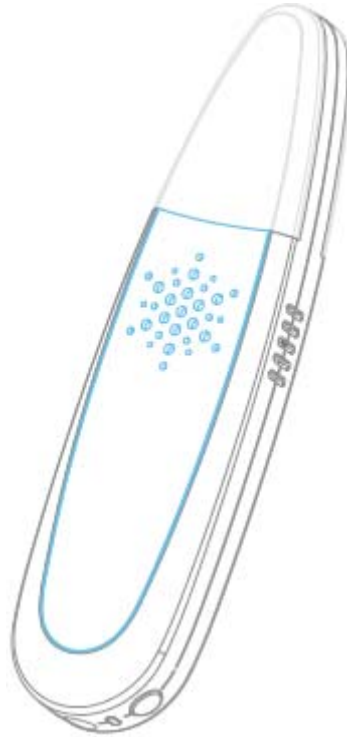


Thanks for your choice of this product.
For suitable and safe use, please understand this manual well to reach the optimal effect in operation.
When this product is transferred to other people, please transfer this service manual with the product and keep it well.

Safety Instructions




Please understand the safety instructions well before operation and use the beauty device according to the operation procedure.

To specify hazard and damage in the following instructions, the details—examples are particularly given so as to help your understanding. The important safety instructions are provided hereof and please follow them in the operation.

The safety instructions are classified into three kinds

 Danger	The misoperation may result in death or severe injuries
 Warning	The misoperation may result in death or severe injuries
 Attentions	The misoperation may result in personal injury and material damages (The material damage means the huge damages such as house, family property and livestock)

Meanings of icons

 Danger	This icon indicates forbidden items
 Warning	This icon indicates items not allowed
 Attentions	This icon indicates the items that the attentions must be paid to

 **Danger**



It is not suitable for the following situations:

Forbid The person who is in driving or in the operation of various electronic equipments

 **Warning**

It is not suitable for the following persons:

1. The person in pregnancy period, menstrual period, growth period and adolescence or with allergic dermatitis, skin sensitivity and allergic constitution.
2. The person who suffers from inflammation and skin fester due to excessive exposure and skin burns.
3. The person who suffers from skin inflammation, skin itching, senile plaques or physical stimulation and other spots due to cosmetics.
4. The person who suffers from viral pigmented naevus and black skin disease.
5. The person with hypertension, weak skin feeling and weak body.
6. The person who uses cholesterol for a long time and suffers from hepatic dysfunction and other disorders.


For
bid

It is not suitable for the following situations or parts:

1. It is not suitable for the patient with acute diseases, tuberculous effusion, febrile disease, infection disease, malignant swelling cut, heart disease or hypertension.
2. It is not suitable for such parts as chest, mouth, private parts and mucosa.
3. It is suitable for the person with physical discomfort, over fatigue, drinking and tired status due to medicine.
4. Please ask a favor from the doctor if you can not make a judgment.


For
bid

The operator must follow the requirements as below:

1. AC adapter must be equipped with AC 100~240V.
2. It is forbidden to use loose socket, extended wire or damaged power plug.
3. Before the adapter is disconnected, the operator must hold the front part of plug but never unplug it with wet hand. In case that abnormal phenomenon such as smoking, abnormal sound and foul smell occurs, the device can not be used. In addition, AC adapter must be disconnected immediately.
4. It is not allowed to use this device in the place nearby combustibles, explosives, fire source or hot appliances, in hot weather or in other places at high temperature.
5. The device is not suitable for the child or the person with mobility, unserved person or those who can not express their thought.
6. It is not allowed to result in failure due to fall or crashing (e.g. falling on the ground).
7. The product shall be kept away from water so it can not be used in bathroom or in bathing or cleaned by hot water.
8. The operator shall not dismount, repair or transform the device arbitrarily. Please call sales store or after-sales service center.


Note

 **Attentions**

The following items must be followed:

1. The device can not be used after exposure. The device shall be used according to the skin status. If there is skin discomfort when the device is used or after the device is used, please stop using the device and see the doctor.



Warni
ng

2. The device can not be operated when the air hole is blocked.

3. Keep the device away from water, impurities and other fluids.

4. The device is not suitable for the person who lies on the bed.

5. It is not allowed to start up the device immediately after shutdown. When automatic running is finished for once or the effect of cold board is not desirable, please wait for cooling or keep the device idle for 1 hour before use.

It is not suitable for the following items:

1. Please use the device according to skin status and stop using the device when the skin suffers from measles, inflammation, itching and other abnormalities.



Warni
ng

2. The device can not be used for other purposes beyond beauty.


3. The device can not be used after falling, crashing or damage.

4. The device can not be used for business or excessive use in continuity.


Contents

Safety Instructions	1
Name of Various Parts	5
Introduction for Use	6
Display Method of Indicator Light	6
Operation Instructions.....	7
Charging Method	7
Reset Method	8
Operation and Curriculum	9
Operation Method	9
Caring Curriculum	11
Operation Completed	14
Maintenance Method	14
Preservation Method	14
Preservation and Inspection	14
Troubleshooting.	14
Doubt Phenomenon.....	14
Common Problems.....	15
After-sales Service	16
Product Specification	16
Warranty Rules.....	16
Guarantee and Declaration.....	18
Guarantee Card.....	18

NAME OF VARIOUS PARTS
PRODUCT APPEARANCES AND NAME OF PARTS

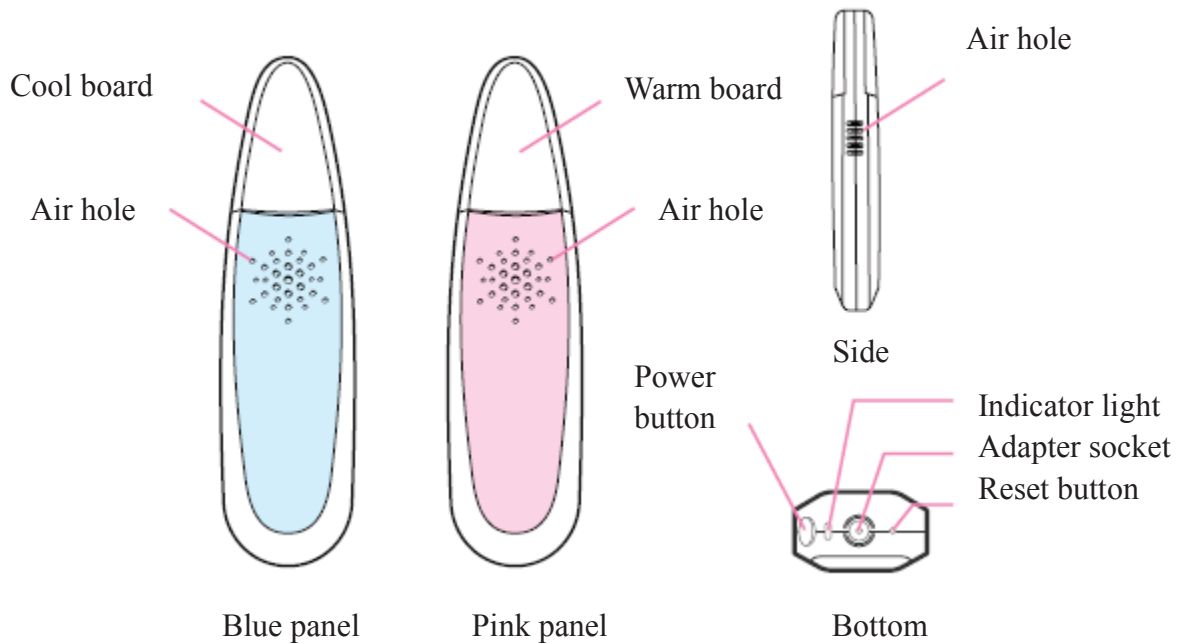

AC ADAPTER

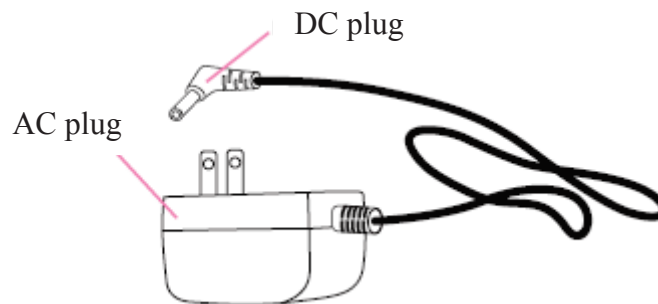

Features of the device include delicate design and simple operation, one button start and stop humanization design shall be adopted.

The operation is simple (Power ON--working mode (10 minutes)--Colding (2 minutes)--Power OFF) shutdown status, keep pressing power button, blue light flickers (i.e. the device starts), when blue light's flicker becomes constant on, the beauty device can be used! After operations keep for 12 minutes, beauty device has auto-timing shutdown function, thus it can shut down automatically. Manual shutdown can also be used, what you need to do is keep pressing power button until the blue light is off (i.e. the device shuts down).

Tips: If you continue to use, please use it again after it is colded down (The interval is recommended to exceed 1 hour), the efficacy shall be better!


INTRODUCTION FOR USE

DISPLAY METHOD OF INDICATOR LIGHT

This machine adopt LED light's flicker (on or off) direction indicate present state of work or charge.

A Indication of light:

Indicator light

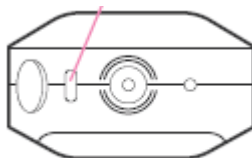

B Display in actions

Blue light is on: Power is on

Press the power button for more than 1 second, the power is on.

Blue light: Enable use state

Indicator light is in flicker of blue light becoming constant on of blue light, that indicate beauty device's temperature is steady, so it can be start used.

Off: Power is off

Press the power button for a while in operation, which can manual shutdown the machine.

C Display in charge

Red light flickers: in charge

Red light is constant on: charge is completed

OPERATION INSTRUCTIONS

Read this operation manual carefully before use and operate according to the manual.

Continue operating time

Please use for no more than 10 minutes in each time (after use the machine may become scorching). For continuous use, please wait until it becomes cool (suggest cooling time is 1 hour).

Air hole

Hold the beauty device, never block air hole, never immerse in the water and prevent the impurities.

Buzz

After startup, the beauty device will buzz. It is not has problem but cooling equipment is cooling. In use, buzz may stop or not.

Care board's temperature

After startup, care board's temperature will be steady in about 90 seconds. Indicator light become blue light's constant on, it can be start used.

※ Room temperature is 25℃

※ Care board's temperature's steady time depends on circumstances.

Care board's operation

After startup, cool board and warm board will work at the same time. The pink side is warm board and blue side is cool board.

You cannot let the machine work single operation, please realize.

CHARGING METHOD


Care board charges by rechargeable battery or adapter. If want to use inner battery, please first make sure inner battery is charged.

Beauty device is lithium battery. Even electric quantity is run out; recharge does not reduce its performance. The first time of using beauty device or long time not use, please charge first and make sure indicator light is flicker of red light become red light's constant on then use it.

CHARGE TIME

Charging time	Operation time
3 hours	60 minutes

- ※ Ambient temperature is 25°C
 - ※ Operation time is time of duration immediately use after charge.
- Adapter DC plug connects to beauty device.

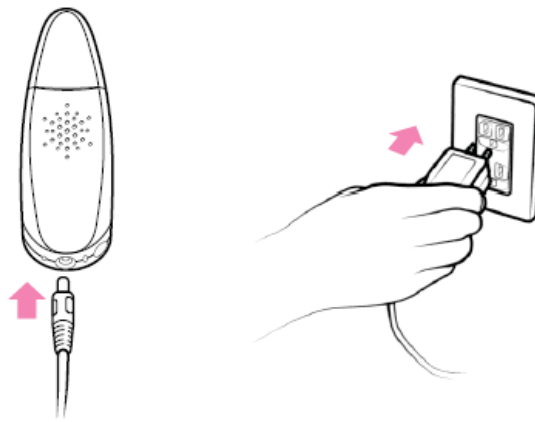

B. Insert the AC plug of the adaptor into the socket.


- When the indicator light flickers red, the beauty device begins to charge.
- After the charging is done, the indicator red light will be on constantly.
- The charging time will be different due to dump energy and ambient temperature difference.
- Press the switch button during charging, and the beauty device will startup, for which, the indicator light will turn to green and then stop charging. The battery can not be charged during operation of the beauty device.

C. After the charging is done, please pull out the AC plug of the adaptor and unplug the AC plug from the machine.

RESET METHOD

When the battery is fully charged, the indicator light won't be on. Please press the "Reset" button and try again.

Reset button. Use a pin to press the button, as shown in the following figure:

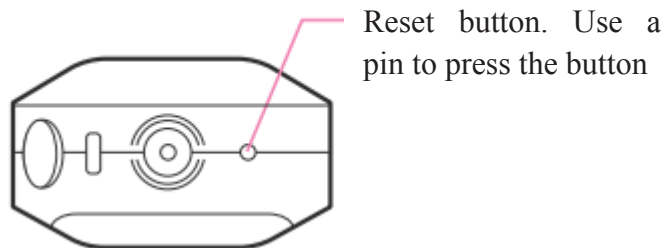

OPERATION & CURRICULUM

OPERATION METHOD


The method to hold beauty device: To use the beauty device as shown in the left figure. Press the press button and the cool board and warm board will begin to work simultaneously. If you want to switch the care mode, please turn the beauty device to the other side.

- Do not block the air hold when holding the beauty device.
- Please touch the cool board and warm board with your finger to feel the temperature before operation.

1. Please press the switch button. The device will startup after pressing the button for 1 second or longer and the indicator light will flicker with green.



Press and tamp down the switch button to startup the device

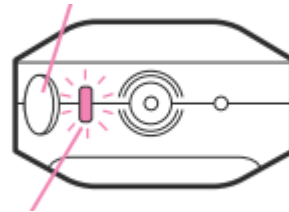


The indicator light will flicker with green


2. Wait the indicator light to become constant green

After 90 seconds. The temperature of the care board will become stable. And the indicator light will become constant green.

- Under a room temperature of 25°C
- The time needed for the care board to become stable depends on factors, such as the ambient temperature, care board temperature before use and the dump energy.

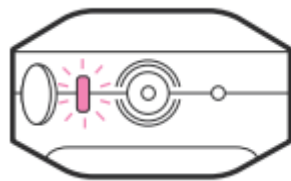

Flickering green

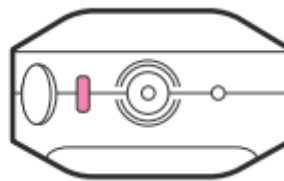

Stable green

3. Conduct cool care with cool board or conduct warm care with warm board.

- The red part refers to the warm board and the blue part refers to the cool board.
- Press or pat the skin gently with the care board. Do not slide it on your skin.

Suggested use: use it for 5 seconds at the same place and turn to another one. 10 minutes per time and at most 5 times a day. Please make proper adjustment according to your own situation.

4. After care

Press the switch button for a while to shutdown the device or it will shutdown automatically after 12 minutes' use.


A. If you want to end the care in the middle of operation, please press the switch button for a while.

B. The warm and cool board has a function of auto shutdown, namely, to enter into stand-by mode, 10 minutes after start. And then, start cooling for about 2 minutes. After the indicator light flickers for about 2 minutes and the beauty device will shutdown automatically.

Shutdown during care, and the about-2-minute cooling function won't start. If you prefer continuous use, please wait for the cooling function to start.



Press and tamp down the switch button to shutdown the device

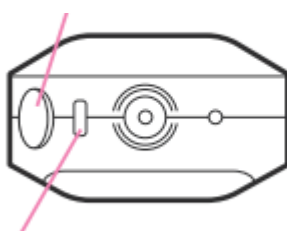

The indicator light will turn off



CAUTION

Tips:

A. After 12 minutes use (include 2 minutes' cooling), the beauty device will shutdown automatically, and the cool board will turn warm slightly. If you want to use it continuously, we suggest you to wait for one hour for an optimum efficiency.

B. Don't block the air hole or ignore the beauty device after startup.


C. Don't use the beauty device with vaporizer or steam vaporizer at the same time. Don't spray water on beauty device or immerse it in water to prevent it from damage.

D. Temperature of the care board may be affected by the ambient temperature. And the warm board may be too hot or too cool, so please feel it with your finger.

CARING CURRICULUM

FIRMING AND SKIN CARE PROCEDURE (3 MINUTES)

Care with warm and cool board in turn to firm and improve the skin.

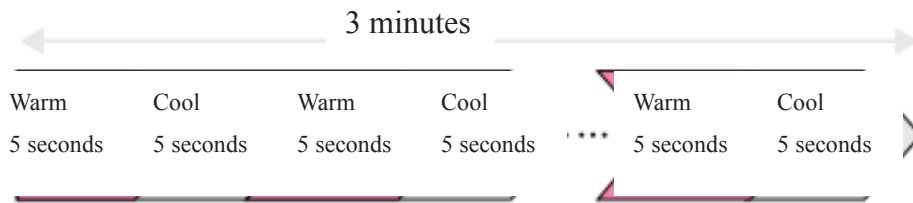

1. To care skin with the warm board for 30 seconds

5 seconds for each part and then turn to the next.

2. Reverse the device and care the same part with the cool board.

2 seconds for each part and then turn to the next.

3. Repeat step 1 and 2 for 3 minutes

Care the whole face. 4-5 seconds for each of the part, such as right cheek, left cheek, wing of nose, jaw line and neck contour and 3 minutes for all.

CARE METHOD


A. To care the cheek, the part below eyes, sides of ears, and edge of lips for a total of 30 seconds. And then, reverse the device to care the same part with the cool board for a total of 15 seconds. Care the left cheek in the same manner after finish the right cheek care.

B. To care the right, middle and left of the forehead with warm board for a total of 30 seconds. And then, reverse the device to care the same part for a total of 15 seconds.

C. To care the wing and bridge of the nose with warm board for a total of 30 seconds and then, care the same part with cool for a total of 15 seconds.

EYE CARE PROCEDURE (1 MINUTE)

Care with warm and cool board in turn to relieve the feeling of strain on your eyes.

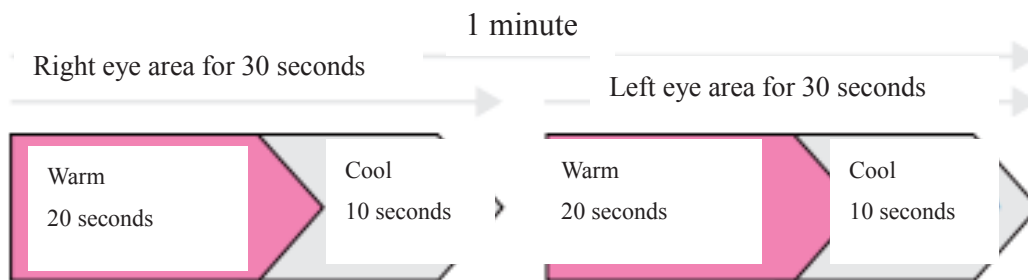

- ※ Care with the warm board from medial canthus to external canthus for 20 seconds
Care for 3 seconds at one part and turn to another one.
- ※ Reverse the beauty device to care with cool board from the part below the eye to area around the eye for 10 seconds.
Care for 2 seconds at one part and turn to another one.
- ※ Care the left eye area in the same manner after finish the right eye area care.
Care for the eye area needs 30 seconds, with a total of 1 minute.


CANTHUS CARE

- A. Care the medial canthus with warm board and then, turn to the external canthus, with a total of 20 seconds.
- B. Reverse the beauty device to care with cool board for 10 seconds.
- C. Use the beauty device to pat the medial canthus gently if you want to care the fatigue eyes. Do not overexert.

BEFORE MAKE UP CARE (3 MINUTES)

Please care with cool board before make up.

Care with cool board for 3 minutes after cleansing, toner, cosmetic emulsion, etc, and then you can make up. As shown in the following figure:

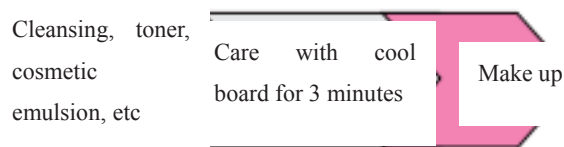

1. Clean your face and do your daily skin care procedure.
Clean your face and apply toner or emulsion in daily manner.
2. Care the face with cool board for 3 minutes.
 - A. Care eye areas, lip area, jaw line, wing of nose and neck contour etc, and each part for 3-5 seconds. And then, turn to another part.
 - B. Press the care board on skin instead of sliding it. And then, uplift it and turn to another part.
 - C. Please refer to the instructions mentioned above for the use of care board.

DEEP CLEAN (3 MINUTES)

Use the warm board to care where there are clogged pores and where there is too much oil. And then, uplift the device and turn to another part. Clean the dirt in deep pores with deep clean and care with warm board for 3 minutes. Clean your face.



1 Use warm board to care clogged pores or parts where oils are too much (3 minutes).

A Care such parts as eye areas, mouth areas, jaw line, wing of nose and neck contour. Every part shall be cared for 3~5 seconds and then transferred to the next part.

B Please refer to aforesaid parts for care board application.

2 Clean the face

Foaming cleanser, other our cleansing massage device and cleansing massage cream you usually used to remove dirt and sebum.

PATTING CARE PROCEDURES (3 MINUTES)

Besides pressing skins, patting method can also be used to care.

The following page describes patting care method in detail. You can also convert previous pages of pressing skin care procedures to patting care method. For patting care, it is appropriate to pat gently, it is not allowed to use excessive force and it is also not applicable for such parts as forehead or nasal bone.

1 Warm board patting care shall be applied on the whole face.

Care board shall be moved simultaneously while patting massage is applied, please refer to following “patting care method”.

2 Cool board patting care shall be applied on the whole face.

Care board shall be moved simultaneously while patting massage is applied, please refer to following “patting care method”.

Patting care method:



A Patting care shall be applied gently along with mouth areas towards the ear direction and then turned to the front of the ear gradually.

B Gentle patting shall be applied along with the nose and then turned to the ear direction gradually.

C Pat the temple gently.


PRODUCTS USED COOPERATIVELY FOR SKIN CARE

It (cleansing or painting skin care products, including beauty cream, toner, beauty fluid, emulsion, vanishing cream, primer potion and BB cream) is optional according to individual demands.

OPERATION COMPLETED

MAINTENANCE METHOD

1. Confirm whether power supply is turned off.
2. Dirt on beauty board is cleaned by wet towel after wringing.

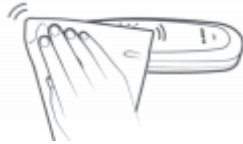

DO NOT use alkaline lotion and alcohol, thinner, gasoline and light removing liquid for using.



- 3 Dirt except on beauty board is wiped by dry cloth.



- 4 Please place it in a place with good ventilation and then preserve it after drying.

PRESERVATION METHOD

When the product is not used, pull out AC adapter from the socket, cut off power supply and preserve it.

DO NOT put it in high temperature and humid environment such as fire source, hot place and hot bath. When the product is not used, please put it in a clean place with good ventilation for preservation.

PRESERVATION & INSPECTION

On condition that this machine is not used for a long time and now it is used again, please firstly confirm the product works in normal condition and then use it.

When abnormality occurs to machine operation, please firstly refer to the following Doubt phenomenon before rework.

If abnormal phenomenon cannot be solved or failure and damage are confirmed, please consult the seller from whom you buy this machine or consult after-sales service center of our Company.

TROUBLESHOOTING

When you feel abnormality when you use it, please confirm whether there is the following phenomenon.

Phenomenon	Possible reasons	Treatment method
Power supply is not electrified	Do not press power switch for a while	Press power switch for a while to turn on power supply
	AC socket is disconnected	Connect AC adapter in a right way according to the requirement
	Rechargeable battery is run out	Please charge it
The product can only be used for several minutes after charging	Charging is insufficient	Please charge it
	Battery life	Please contact after-sales service center of the company
Indicator red light does not flicker or constant on does not occur	AC adapter is not inserted in the socket	Please insert AC adapter
	Power plug is not inserted in beauty device	Insert plug is inserted in beauty device
	Battery is not charged	Please continue charging
Temperature of cool board rises	Continuous use for more than 10 minutes	Please wait for cooling or use it again after placing for one hour
Temperature of cool board will not fall		
Temperature of warm board will not rise		
Temperature of warm board is higher than that of	Air hole is blocked or placing it in blocking state	Use it again when air hole is not blocked after one hour

COMMON PROBLEMS

Q Can beauty cream and emulsion can be used during operation?

A As beauty cream and others will cause failure when they enter inside of beauty device, thus, it is not allowed to immediately use beauty board when directly applying beauty cream and others on the board or applying them on the skin. Beauty cream and emulsion shall be applied on the skin before operation or after operation completion.

Q Can it be used during having a bath?

A Beauty device does not have waterproof function. It cannot be used in bathroom and the place with water and high temperature layer.

Q Can it be used with portable vaporizer and atomizer?

A Although it can be used before and after operation, it cannot be used with portable vaporizer and atomizer. That is because inside of this machine may accumulate moisture, which will cause failure.

AFTER-SALES SERVICE

Rechargeable battery is consumable and it does not belong to warranty object.

If the battery is used for a very short period after charging, it is estimated to reach service life. The target charging time is 300 times.

When you want to change rechargeable battery, please consult sales store or after-sales service point.

Product specification

1	Power supply	DC 5V 1.8A
2	Power	About 12.5W
3	Battery for use	Lithium battery 3.7V/1200 MAH
4	Charging time	About 3 hours
5	Product weight	About 175G
6	Overall dimension	Length about 174*40*25MM

WARRANTY RULES

When failure occurs during normal use, we can implement one year warranty. Consumables are excluded. You do not need to pay any maintenance cost within warranty period. However, you need to pay transport charge.

When failure occurs beyond warranty period, you need to pay for repair. You also need to pay for repair within warranty period under some conditions. So please carefully read instruction for free repair rules.

In case of failure, please consult the seller or after-sales service point.

GUARANTEE

1 When failure occurs to attentions such as label attached on beauty device in normal use, you can receive free repair in point of sales according to service manual.

2 When failure occurs within free warranty scope during warranty period, please receive free repair in point of sales with this machine and the guarantee.

3 Please tell us your detailed mailing address for returning.

4 If free product in record cannot be repaired, please directly contact our company.

5 You have to pay for repair during warranty period in the following occasions.

A The failure and damage caused by wrong using method, disassembly by oneself, wrong repair and refit

B The failure and damage caused by falling after purchase and falling down on the ground

C The failure and damage caused by fire, earthquake, flood and other natural disaster

D The failure and damage caused by connecting another power supply or charger

E The failure and damage of components such as lithium battery and consumables

F Except ordinary family, the failure and damage caused by long time operation

G The failure and damage caused by transportation, handling and falling after purchase

H The occasion which is not listed in the guarantee clearly

I The date, customer name and point of sales are not filled in the guarantee or there is no seal.

J The guarantee is valid in China.

6 The guarantee is applied to free repair of beauty device with clear time and meeting conditions. Other conditions are restricted to some extent. If you need repair beyond warranty period, please consult the seller. For instance, battery consumable needs to be changed.

DECLARATION

For better service, the content of service manual will be changed with product upgrading in time without further notice.

If practical operation of the product is not consistent with service manual, please visit company website and inquire the latest product information.

Our company will not assume any legal liability for the loss caused by using this machine or requirements of various severe conditions. All power of interpretation is reserved by our company.

GUARANTEE CARD

NO			
Goods		Model	
Product name			
Customer name		Tel	
Address			
Purchase date			
Free repair	Object	Complete machine	Exclusive of consumables
	Period	One year	Calculated from the purchase date
Sales store name			
Repair record I			
Repair record II			
Repair record III			
Warranty rules	Please see warranty notes.		